

STAIRCASE MEASURING GUIDE

for the UP Curved Stair Lift

In order to ensure that your AmeriGlide UP Curved Stair Lift installation goes smoothly, it is imperative that proper and accurate measurements are provided. If you have any questions, or would like some assistance with measuring your staircase, please give us a call at **1-800-790-1635**.

1. Photograph your Complete Staircase

Please take photos of your entire staircase from the bottom, middle and top. Please be sure to take as many photos as possible and include any obstructions such as doorway and nearby walls. Please email or text the photos to your AmeriGlide Sales Representative.

2. Staircase Width

inches

Measure the width of your staircase at the most narrow point, then enter the measurement.

3. Rise of One Step*

inches

Measure the height from the top of a step to the very next step down.

4. Tread Depth of One Step*

inches

Measure from the front of one step to the front of the kickplate.

5. Bottom Clearance

inches

If there is a wall or a door at the bottom of the steps measure from the bottom step to the wall.

6. Mid-Landing Width

inches

Measure the width of your staircase landing and enter the measurement.

7. Mid-Landing Depth


inches

Measure the depth of your staircase and enter the measurement.

8. Top Clearance (of final flight)

inches

Measure from the top step to the nearest door, wall, or other obstruction at the top of the stairs.

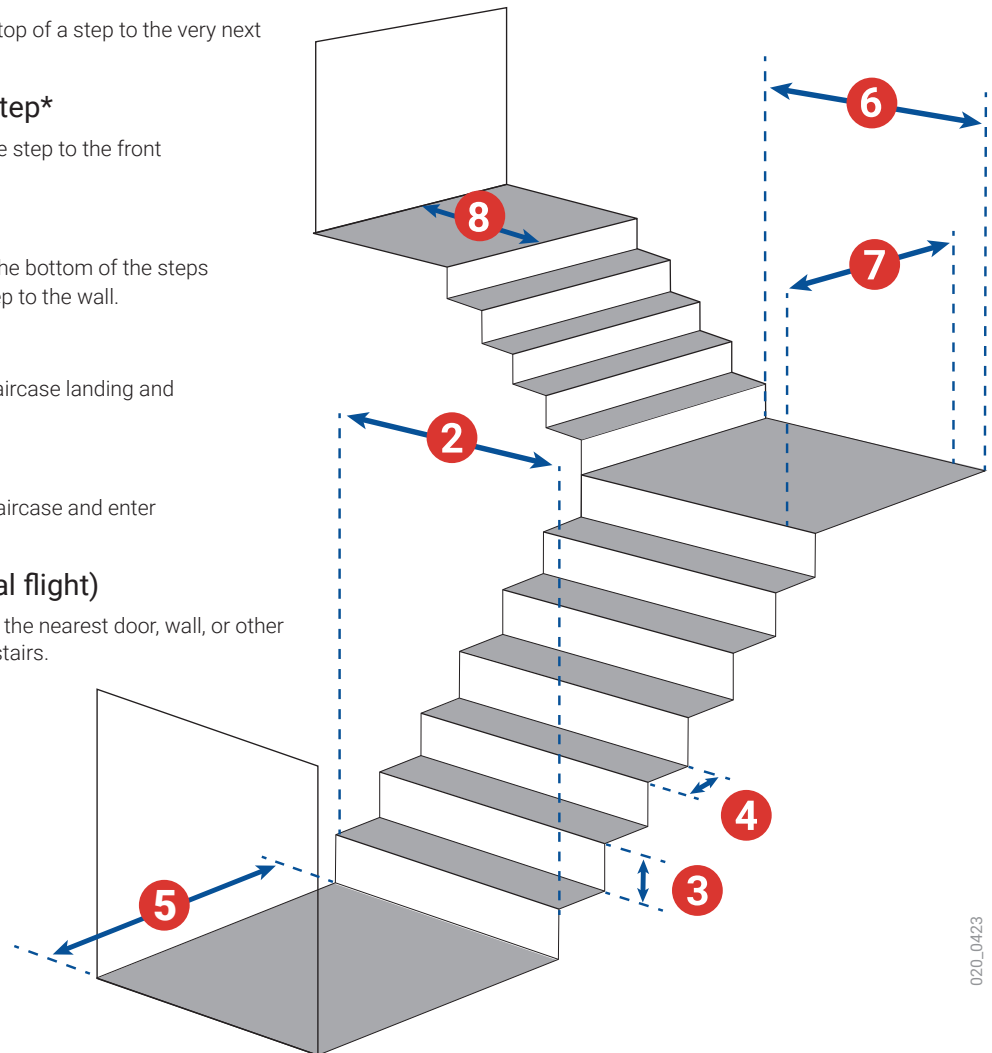


***NOTE:** If your staircase includes different shaped steps, we will need to obtain the measurements for each step type.

Need Assistance?
1-800-790-1635

In order to properly configure your stair lift for your staircase, please provide the information below. If at any point you have any questions, or would like some assistance with measuring your staircase, please give us a call at **1-800-790-1635**.

9. How will your Stair Lift be Oriented on the staircase?

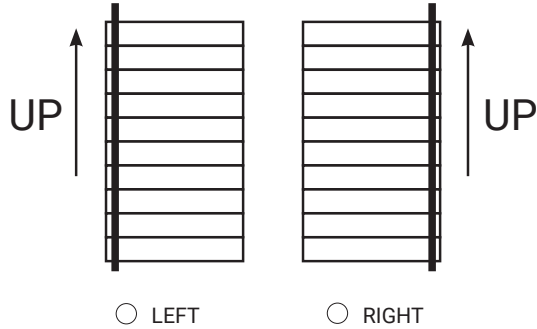

11. Floor Construction:

- ☐ WOOD ☐ CONCRETE ☐ TILE
☐ MARBLE ☐ GRANITE ☐ STEEL
☐ OTHER _____


12. Power Supply (socket) Location:

- ☐ TOP ☐ BOTTOM ☐ MID LANDING


10. Total Number of Stops _____

13. Total Number of Steps _____


14. How would you like your Stair Lift to park at the **BOTTOM** of the staircase?

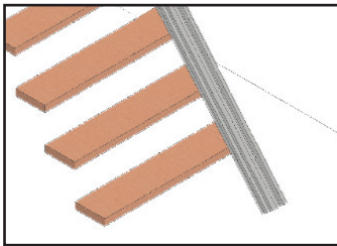


☐ STRAIGHT

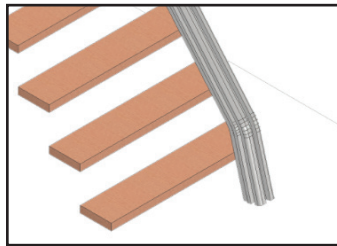


☐ DROP

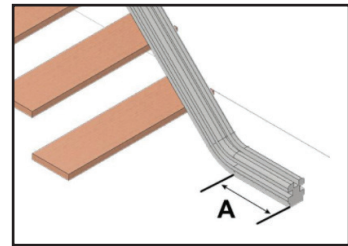


☐ HORIZONTAL

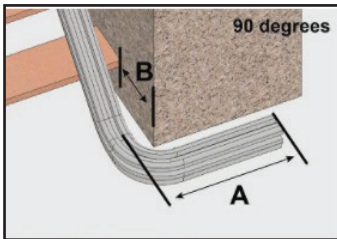

☐ 90° Parking

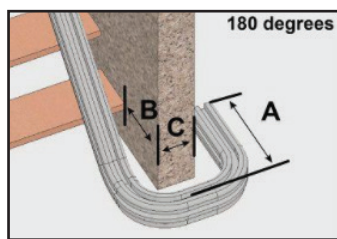


☐ 180° Parking

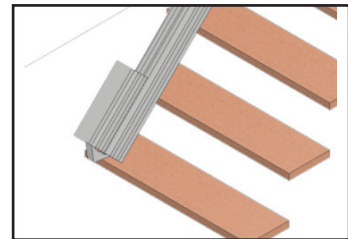


☐ HINGE RAIL

Need Assistance?
1-800-790-1635


15. How would you like your Stair Lift to park at the **TOP** of the staircase?

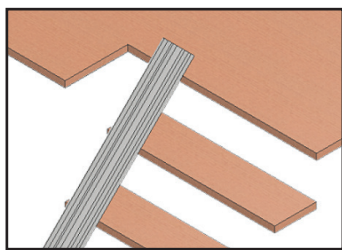


☐ STRAIGHT

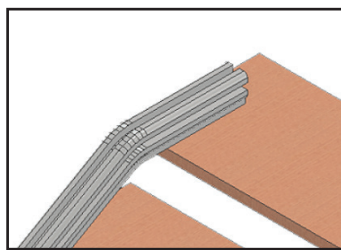

☐ HORIZONTAL

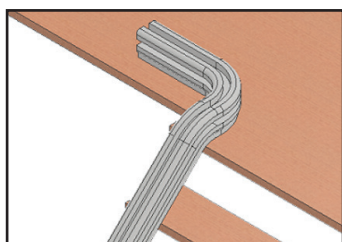

☐ 90° Parking

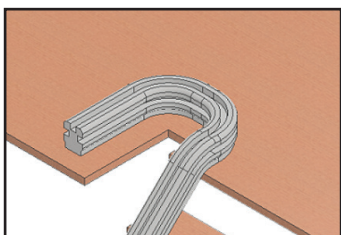

☐ 180° Parking

Continue to next page.

Need Assistance?
1-800-790-1635

**AmeriGlide**
ACCESSIBILITY SOLUTIONS